





We maximize customer values with E-CUBE principles



The ALPINION E-CUBE Series continues to evolve through unique and reliable technological breakthroughs.

The E-CUBE 9 Series integrates ALPINION's core imaging technologies to provide uniform and fundamentally excellent image quality across the whole product lifetime.

E-CUBE 9 Professional



E-CUBE 9 DIAMOND



The outstanding quality of E-CUBE 9
Series results from *a refined acoustic technology* with great delicacy that
manages and handles ultrasound signals
perfectly.

The ALPINION's unique and flexible imaging platform, FleXcanTM
Architecture, enhances system stability and helps to maintain system performance to the latest software version. The reliable imaging platform guarantees uniform image quality.

The unique imaging platform and proprietary acoustic technology of E-CUBE 9 Series delivers *uniform image quality and crystal clear image performance* to clinicians.

The evolution of ALPINION's core imaging technology

ALPINION's innovative imaging platform (FleXcan[™] Architecture) and sophisticated image processing technology set (Optimal Imaging Suite[™]) provide a consistent and high-quality image performance that is ideally suited to a dynamically changing work environment.

FleXcan[™] Architecture

The ALPINION's unique ultrasound imaging platform (FleXcan[™] Architecture) is a software driven structure that ensures stable imaging performance, easy software upgrades and uniform image quality throughout the whole product lifetime.

Optimal Imaging Suite[™]

The ALPINION's integrated post-processing technology set (Optimal Imaging Suite™), embedded in the E-CUBE 9 Series, creates optimized images for each application by effectively decreasing artifacts and enhancing the borderline of each organ.

Crystal Signature™

The ALPINION's unique acoustic engineering produces a superior beam-field pattern that provides a uniform image performance. ALPINION's specialized single crystal technology produces a characteristically sharp and deep elevation beam focus.

The innovative single crystal technology is characterized by a higher energy conversion efficiency than conventional pizeo-ceramic materials that provides greater uniformity and sensitivity.


The focused beam pattern, produced by ALPINION's single crystal technology (Crystal Signature[™]), results in high-resolution image, wide-bandwidth images, and deep penetration.

[Speckle Reduction Imaging technology Spatial Compounding technology Beam Field (normalized at each depth) Competitor ALPINION


SensitiView[™]

The ALPINION's unique and proprietary acoustic technology,

SensitiViewTM, encompasses the CSATM (Clear Signal Amplifier), low loss cable and matching technology with system. The CSATM generates an enriched purified signal with active electronics from high-quality piezo-electric materials. The high-performance low loss cable transfers high-sensitive signals for optimum impedance matching with the system. The results is significantly enhanced sensitivity and high-resolution images.

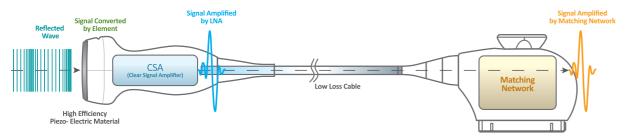

Innovation is born when design cares daily practice

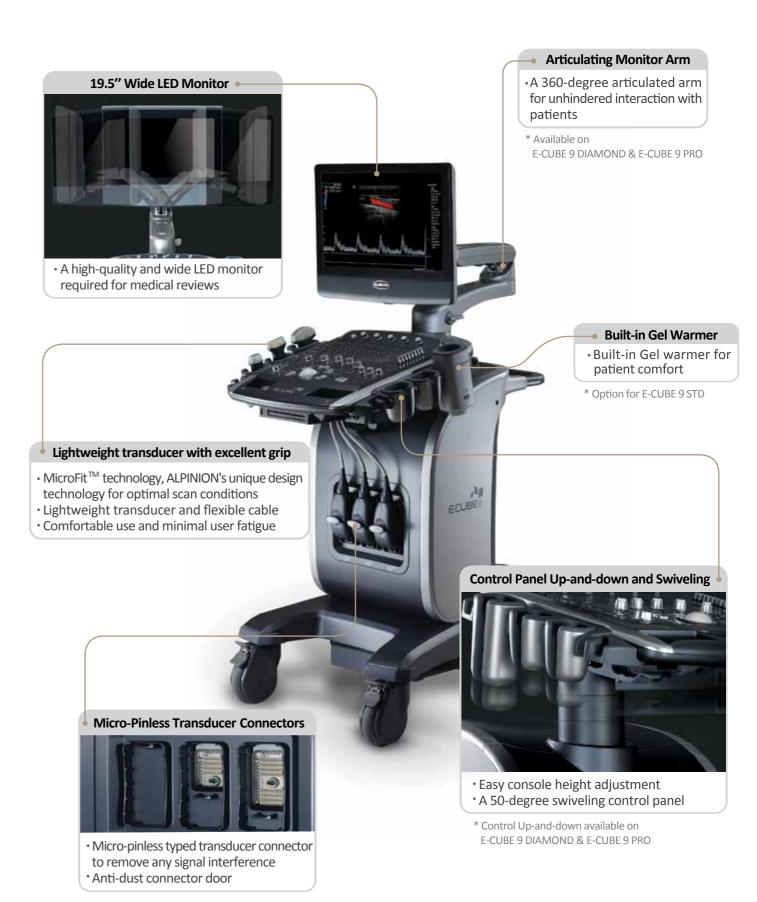

Efficient diagnostic flow

The efficient workflow adds to the performance perfection of E-CUBE 9 Series.

The E-CUBE 9 Series enhances user value in each phase of daily practice, from patient registration to clinical review, as well as diagnostic confirmation. The E-CUBE 9 Series fine-tunes its imaging technology, dedicated application software and built-in features for an optimized diagnostic flow.

Clinician's value:

• Minimizes the time and effort to register a patient.

Solution • Quick ID

Simplifies patient registration process by basic information only. (Patient ID, name, year of birth and gender).

Clinician's value:

- Application-specific software to maximize diagnostic confidence
- Easy operation through intuitive user interface

Solution • Auto IMT

Carotid scanning observes early stage arteriosclerosis.

The ALPINION's Auto IMT software enables accurate, automatic and precise measurements of the thickness of the upper and lower intimamedia through a single line near the intima.

Live HQ™

The combination of qualified 3D/4D transducers with strong performance reliability of geometrical accuracy enables ALPINION's Live HQ^{TM} software to render a realistic fetal shape. The innovative illumination effect enhances the diagnostic confirmation of a morphological examination.

CUBE Strain™

The ALPINION's strain imaging enables the quantitative and objective evaluation of myocardial contractions. CUBE Strain™ detects velocity gradients and helps clinicians analyze regional myocardial function.

Needle Vision™ Plus

The ALPINION's sophisticated needle visualization technology, Needle Vision™ Plus dramatically enhances visibility of needle tip by increasing perpendicular reflection with Beam steering technology.





• Minimizes the time and effort to register a patient.

Solution • Power User Preset

Quick mode transition based on pre-specified user presets that include probe selection, assigned specific function and image parameters.

Clinician's value:

• Image optimization with minimal manipulation.

Solution • Xpeed™

Xpeed[™] optimizes images in B-mode, color flow and Doppler mode with one button operation.

FullSRI™

The ALPINION's unique speckle reduction filtering technology (FullSRI[™]) removes artifacts and enhances tissue edges with a high-resolution contrast.

SCI (Spatial Compounding Imaging) and FCI (Frequency Compounding Imaging)

The SCI and FCI provides higher contrast, less artifacts and more uniform images (from near field to far field) by steering multiple beam angles and combining ultrasound data from different frequencies.

Defining pure and high resolution images

Optimized for abdominal scans

- Excellent performance in B/C/Doppler mode
- Deeper penetration with high-sensitivity image from Crystal Signature[™] technology
- A wide selection of transducers that include various high-density convex transducers




- Superior B/C/Doppler performance
- Clear needle detection imaging
- A wide range of linear transducers



The ALPINION's excellent linear transducer performance is a signature feature of E-CUBE 9 Series. Its qualified image performance provides clear needle detection images for accurate FNA (Fine Needle Aspiration) and needle injection.

The E-CUBE 9 Series offers the widest range of linear transducers, from the hockey stick to the wide footprint linear transducer. This optimal transducer portfolio assists in the clinician's diagnostic specialty.

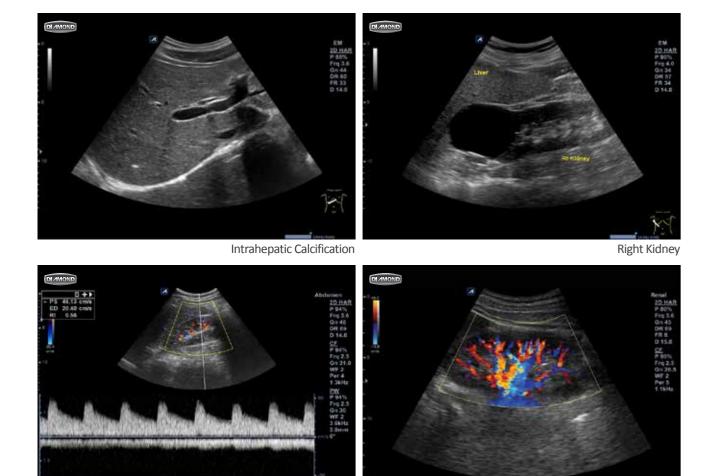

Renal Spectrum


Renal Color Flow

The E-CUBE 9 Series features advanced single crystal transducers for deep penetration, wider bandwidth

Found in high-end systems, the high-density convex transducer of the E-CUBE 9 Series heightens diagnostic

and high-sensitivity images. The results are exceptional even from high-volume patients.

confidence and offers uniform high-resolution images in every abdominal examination.

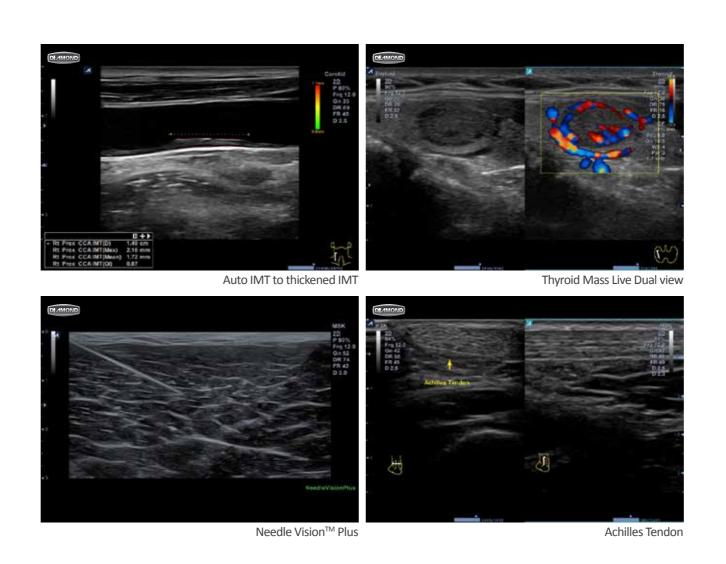

Optimized solutions for precise women's health exams


- Intuitive and lucrative 3D/4D software (Volume Master™)
- Superb and realistic fetus images rendered by Live HQ[™]
- Full transducer portfolio for women's health, from breasts to 3D/4D



The E-CUBE 9 Series provides optimized clinical solutions for all types of women's health exams through outstanding image performance and intuitive workflow.

The ALPINION's technological breakthroughs create the world-first single crystal volume convex transducer and the high performance 60mm wide linear transducer for breast exams.

Volume Master[™], ALPINION's intuitive 3D/4D software, offers fast rendering speed and outstanding rendering quality. The innovative rendering technology, Live HQ[™], provides realistic fetus images and fosters maternal-fetal bonding.

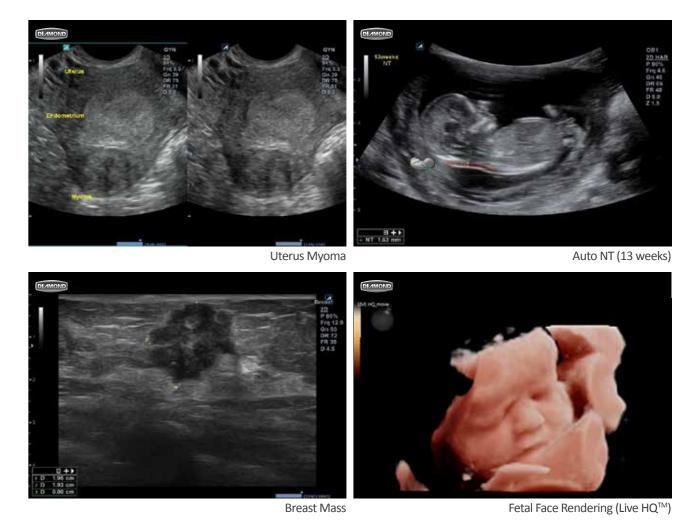

Customized package for cardio-vascular exams

- Application dedicated transducer sets, from pediatric to adult cardiac exams
- Accurate measurement packages and report function
- Comprehensive CV tools and software to support diagnostic decisions



Integrating ALPINION's core imaging technology, a wide range of CV specialized transducer sets, accurate measurement, comprehensive CV tools and structured connectivity, the E-CUBE 9 Series offers clinicians a total diagnostic cardio vascular examination solution.


The E-CUBE 9 Series provides a wide selection of transducer sets, from single crystal adult and pediatric cardiac transducers to the peripheral deep vascular-specialized transducer.

The full CV measurement package and CV-specialized tools (ECG, CUBE Strain™, Stress Echo and Auto IMT) support more precise cardio-vascular examination. A streamlined workflow with the seamless connectivity of DICOM 3.0, the modality worklist and the DICOM Cardiac SR guarantee clinicians reliable diagnostic efficiency in routine cardiac examination.



Tissue Doppler Imaging

Tricuspid Regurgitation

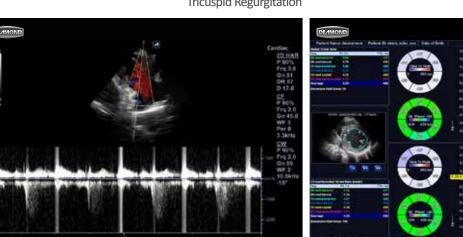

Tricuspid Regurgitation CW

CUBE Strain^{TI}

Premium class transducers

The ALPINION's uniform and fundamentally excellent image performance is the result of an unique and proprietary acoustic technology. The E-CUBE 9 Series's premium class transducers provide clinicians with high definition image, excellent grip, easy contact and reliable performance.



Developed and manufactured by ALPINION







* A biopsy kit is available.

Convex



High dengity Single Crystal

- OB, GYN, Fetal Echo, Abdominal, Renal, **Emergency Medicine**
- * Available for E-CUBE 9 DIAMOND only



Wide-angle high density Single Crystal

- Abdominal, Renal, OB, GYN, Fetal Echo, Emergency Medicine
- * Available for E-CUBE 9 DIAMOND only



High density Single Crystal

• Abdominal, Renal, OB, GYN, Fetal Echo, **Emergency Medicine**



for unrivalled acoustic performance

SC1-6 * Single Crystal

• Abdominal, OB, GYN, Fetal Echo, **Emergency Medicine**



• Abdominal, OB, GYN, Fetal Echo, **Emergency Medicine**



C5-8N / C5-8 Micro convex

• Pediatric, Cardiac, Vascular, Emergency Medicine

Volume Convex



High dengity Single Crystal

OB. GYN. Fetal Echo



VC1-6

OB, GYN, Fetal Echo



L8-17H High density

• Carotid, Peripheral Vascular, Thyroid, Testicle, MSK, Superficial, Breast, Emergency Medicine



• Carotid, Peripheral Vascular, Thyroid, Testicle, MSK, Superficial, Breast, Emergency Medicine

L8-17



High density

• Carotid, Peripheral Vascular, Thyroid, Testicle, MSK, Breast, Emergency Medicine



HD Wide linear array (64mm footprint)

Breast, Thyroid, MSK

L3-12 * Linear array

 Carotid, Peripheral Vascular, Thyroid, Testicle, MSK, Breast, Emergency Medicine

Endocavity



Low frequency with the greater depth of penetration

• Peripheral Vascular, MSK, Emergency Medicine



High frequency hockey stick

• MSK(Superficial), Small Parts, Superficial,



103-12 Intraoral

• Salivary gland, Parotid gland, Submaxillary gland, Small Parts, MSK, Peripheral Vascular, Emergency Medicine



Single crystal phased array for adult

 Cardiac, Abdominal, Renal, TCD, Emergency Medicine



Single crystal phased array for pediatrics

• Abdominal, Cardiac, Emergency Medicine



• GYN, OB, Fetal Echo, Emergency Medicine

EV3-10H * High density, wide bandwidth

EC3-10H * High density, wide bandwidth · GYN, OB, Fetal Echo, Urology,

Emergency Medicine



EC3-10 *

• GYN, OB, Fetal Echo, Urology, Emergency Medicine



EV3-10 *

• GYN, OB, Fetal Echo, Emergency Medicine



Pencil typed transducer

Cardiac



Pencil typed transducer

Cardiac



Real Time Customer Service

NETCUBECARE[™]

ALPINION's E-CUBE system can be accessed by Net CUBE CARE™ service anytime for the efficient troubleshooting of system software. Net CUBE CARE™, ALPINION's unique network-based remote service, efficiently upgrades the E-CUBE system based on customer needs by accessing and upgrading software. ALPINION's Net CUBE CARE™ customer service is optional - please contact the ALPINION customer service team at service@alpinion.com for further inquiries.

ALPINION MEDICAL SYSTEMS

Verdi Tower, 72, Digital-ro 26-gil, Guro-gu, Seoul, 08393 Korea

Homepage : www.alpinion.com E-mail : international@alpinion.com

TEL:+82-2-3282-0900 FAX:+82-2-851-5591

Standalone clinical images may have been cropped to better visualize pathology.

Copyright©2015 ALPINION MEDICAL SYSTEMS Co., Ltd. All rights reserved. Catalogue contents may be partly changed without prior notice to customers due to performance enhancements.